

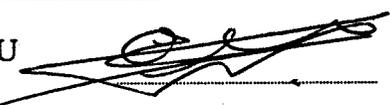
DRAFT SIGN OFF BY THE WORKING GROUP
STEP 2 DOCUMENT

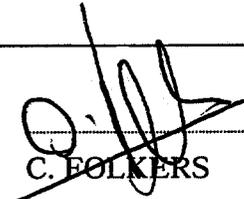
Topic Reference VICH GL _____

Subject Biological Quality Monitoring Working Group
: Residual Formaldehyde Testing

Draft No. 3.2a Dated 07 / 11 / 00 _____

EUROPE

EU 
(NAME) R. DOBBELAER

FEDESA 
(NAME) C. FOLKERS

JAPAN

JMAFF 伊藤 治
(NAME) O. ITOH

JVBA 藤川 勇治
(NAME) Y. FUJIKAWA

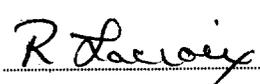
USA

USDA 
(NAME) P. F. ROSS

AHI 
(NAME) H. DRAAYER

OTHER EXPERTS

AUSTRALIA/NEW ZEALAND 
(NAME) F. BOVER

CANADA 
(NAME) R. LACROIX

VICH, International Cooperation on Harmonization of Technical Requirements for Registration of Veterinary Medicinal Products

WORKING GROUP: Biologicals Quality Monitoring

TOPIC: Test for Residual Formaldehyde

GUIDELINE
FOR TESTING OF
RESIDUAL FORMALDEHYDE
(v3.2a)

TABLE OF CONTENTS

INTRODUCTION

1. Objective of Guideline
2. Scope of the Guideline
3. Background
4. General Principle

PROPOSED GUIDELINE

1. Definition of samples to be used for testing
2. General Test procedure

APPENDIX

1. Formaldehyde level conversion table
2. Current regulatory considerations on product residual formaldehyde

REFERENCES

INTRODUCTION

1) Objective of Guideline

Many inactivated veterinary vaccines, particularly bacterins contain residual levels of formaldehyde. It is important to determine the residual level of formaldehyde to:

- a) help assure product safety,
- b) assure the product will not inactivate other products used in combination,
- c) help assure the product remains active throughout its shelf life, and
- d) assure any clostridial toxoids will be antigenic and safe.

This document provides a guideline for the general requirements for residual formaldehyde testing. The guideline leaves the flexibility for other testing methods based on the specific scientific situations or characteristics of the target material. These variations must be stated in the manufacturers production method and include equivalence data.

2) Scope of Guideline

This guideline applies to final product testing for formaldehyde containing veterinary vaccines.

3) Background

A number of assays are available for the determination of residual free formaldehyde in inactivated vaccines, including acetyl acetone titration, ferric chloride titration and the basic fuchsin test. The ferric chloride method was selected since it has shown to be compatible with products neutralized with sodium bisulphite.¹

Residual formaldehyde will be reported as g/L, a conversion table is provided as Appendix 1. Current regulatory considerations on product residual formaldehyde levels are provided in Appendix 2.

4) General Principle

Total formaldehyde is determined based on the reaction of formaldehyde with Methylbenzothiazolone hydrazone hydrochloride (MBTH). The method involves (a) the combination of MBTH and formaldehyde to give one product, (b) the oxidation of excess MBTH to give another product and (c) the combination of these two to give a blue chromophore which is measured at 628 nm.²

Ferric Chloride Method

1. Reagents

- 1.1. *Ferric chloride-sulphamic acid reagent.* A solution containing 10 g/L of *ferric chloride* and 16 g/L of *sulphamic acid*.
- 1.2. *Methylbenzothiazolone hydrazone hydrochloride reagent.* (MW 233.7). [CAS 149022-15-1]. *3-Methylbenzothiazol-2(3H)one hydrazone hydrochloride monohydrate.* An almost white or yellowish crystalline powder. mp: about 270 °C. A solution containing 0.5 g/L L (WARNING: This solution is not stable and should be prepared fresh daily).
- 1.3. *Suitability for determination of aldehydes.* To 2 ml of *aldehyde-free methanol* add 60 µl of a 1 g/L solution of *propionaldehyde* in *aldehyde-free methanol* and 5 ml of a 4 g/L solution of *methylbenzothiazolone hydrazone hydrochloride*. Mix, allow to stand for 30 min. Prepare a blank omitting the *propionaldehyde* solution. Add 25.0 ml of a 2 g/L solution of *ferric chloride* to the test solution and to the blank, dilute to 100.0 ml with *acetone R* and mix. Measure absorbance of the test solution on a spectrophotometer at 660 nm in a 1 cm cell using the blank as compensation liquid. The absorbance of the test solution must be greater than or equal to 0.62 absorbance units.
- 1.4. *Formaldehyde solution,* containing not less than 34.5 percent w/v and not more than 38.0 percent w/v of formaldehyde (CH₂O)
- 1.5. *isopropyl myristate,* analytical grade
- 1.6. *hydrochloric acid* (1 M), analytical grade
- 1.7. *chloroform,* analytical grade
- 1.8. *sodium chloride* (9 g/L and 100 g/L aqueous solutions), analytical grade
- 1.9. *polysorbate 20,* analytical grade

2. Sample and Standards Preparation

- 2.1 Prepare *formaldehyde* standards of 0.25, 0.50, 1.00 and 2.00 g/L by diluting *formaldehyde solution* (1.3) with water in suitable volumetric flasks.
- 2.2 If vaccine to be examined is an oil emulsion, the emulsion should be broken by a suitable method. The formaldehyde concentration in the aqueous phase should be measured. The following separation techniques have been shown to be appropriate.
 - a) Add 1.00 ml of vaccine to 1.0 ml of *isopropyl myristate* and mix. To the mixture, add 1.3 ml of 1 M *hydrochloric acid*, 2.0 ml of *chloroform* and 2.7 ml of 9 g/L *sodium chloride*. Mix thoroughly. Centrifuge at 15,000 g for 60 min. Transfer the aqueous phase to a 10 ml volumetric flask and dilute to volume with water. Use the diluted aqueous phase for the test for formaldehyde. If this procedure described fails to separate the aqueous phase, add 100 g/L of *polysorbate 20* to the sodium chloride solution and repeat the procedure, but centrifuge at 22,500 g.

- b. Add 1.00 ml of vaccine to 1.0 ml of a 100 g/L solution of sodium chloride and mix. Centrifuge at 1000 g for 15 minutes. Transfer the aqueous phase to a 10 ml volumetric flask and dilute to volume with water. Use the diluted aqueous phase for the test for formaldehyde.
- c. Add 1.00 ml of vaccine to 2.0 ml of a 100 g/L solution of sodium chloride and 3.0 ml of chloroform and mix. Centrifuge at 1000 g for 5 minutes. Transfer the aqueous phase to a 10 ml volumetric flask and dilute to volume with water. Use the diluted aqueous phase for the test for formaldehyde.

NOTE: Volumes used for breaking emulsions are for the purpose of illustration. Volumes may differ subject to proportional adjustment of the volumes of other reagents used in the extraction process.

3.0 Test Method

- 3.1 To 0.50 ml of a 1 in 200 dilution of the vaccine to be examined (if emulsion, use 0.50 ml of a 1 in 20 dilution of the diluted aqueous phase), and to 0.50 ml of 1 in 200 dilution of each of the formaldehyde standards, add 5.0 ml of the methylbenzothiazolone hydrazone hydrochloride reagent. Close the tubes, shake, and allow to stand for 60 min.
- 3.2 Add 1 ml of ferric chloride-sulphamic acid reagent and allow to stand for 15 min.
- 3.3 Measure absorbance of vaccines and standards on a spectrophotometer at 628 nm in a 1 cm cell, using the reagent blank as compensation liquid.

4.0 Calculations and Interpretation

Calculate total formaldehyde concentration (g/L) from the standard curve using linear regression

Appendix 1

Formaldehyde level conversion table

g/L formaldehyde	% w/v formaldehyde	% v/v formaldehyde solution*	ppm formaldehyde
2.0	0.2	.5	2000
0.8	0.08	0.2	800
0.4	0.04	0.1	400
0.5	0.05	0.125	500
0.05	.005	0.0125	50
0.04	.004	0.01	40

* based on 40% formaldehyde solution

Appendix 2

Residual Formaldehyde concentration (g/L)	Safety considerations
2.0	Upper limit, but higher levels may be acceptable on a case by case basis and will be supported by appropriate safety data with product assembled at the approximate maximum residual formaldehyde level, including the evaluation of injection site reactions.
0.5 to 2.0	Target animal safety data may be required, with product assembled at the approximate maximum residual formaldehyde level, including the evaluation of injection site reactions.
< 0.5	No additional safety data required.

References

1. Chandler, M.D. & G.N. Frerichs, Journal of Biological Standardization (1980) 8, 145-149
2. Knight, H, & Tennant R.W.G. Laboratory Practice, (1973) 22, 169-173