

United States Department of Agriculture
Center for Veterinary Biologics
Testing Protocol

SAM 213

Supplemental Assay Method for Potency Testing of
Clostridium botulinum Type C Bacterin-Toxoids

Date: May 17, 2016

Number: SAM 213.03

Supersedes: SAM 213.02, December 28, 2009

Standard Requirement: 9 CFR 113.110

Contact: Janet M. Wilson, (515) 337-7245

Approvals: /s/Larry R. Ludemann Date: 01Jun16
Larry R. Ludemann, Section Leader
Bacteriology

/s/Paul J. Hauer Date: 08Jun16
Paul J. Hauer, Director
Policy, Evaluation, and Licensing
Center for Veterinary Biologics

United States Department of Agriculture
Animal and Plant Health Inspection Service
P. O. Box 844
Ames, IA 50010

Mention of trademark or proprietary product does not constitute a guarantee or warranty of the product by USDA and does not imply its approval to the exclusion of other products that may be suitable.

| |
|--|
| Entered into CVB Quality Management System by: <u>/s/Linda S. Snavelly</u> <u>13Jun16</u> Linda S. Snavelly Date Quality Management Program Assistant |
|--|

UNCONTROLLED COPY

Supplemental Assay Method for Potency Testing of *Clostridium botulinum* Type C Bacterin-Toxoids

Table of Contents

- 1. Introduction**
- 2. Materials**
 - 2.1 Equipment/instrumentation**
 - 2.2 Reagents/supplies**
 - 2.3 Test animals**
- 3. Preparation for the Test**
 - 3.1 Personnel qualifications/training**
 - 3.2 Selection and handling of test animals**
 - 3.3 Preparation of supplies**
 - 3.4 Preparation of reagents**
- 4. Performance of the Test**
 - 4.1 Vaccination**
 - 4.2 Challenge of immunity**
- 5. Interpretation of the Test Results**
- 6. Report of Test Results**
- 7. References**
- 8. Summary of Revisions**

Appendix

UNCONTROLLED COPY

Supplemental Assay Method for Potency Testing of *Clostridium botulinum* Type C Bacterin-Toxoids

1. Introduction

This Supplemental Assay Method (SAM) describes procedures for potency testing biological products containing *Clostridium botulinum* Type C antigen as prescribed in title 9, *Code of Federal Regulations* (9 CFR), part 113.110. Mink are vaccinated subcutaneously then challenged intraperitoneally 21 to 28 days later with a standard dose of *C. botulinum* Type C toxin.

2. Materials

2.1 Equipment/instrumentation

- 2.1.1 Vortex-type mixer
- 2.1.2 Refrigerator, 2°- 7°C
- 2.1.3 Freezer, -70°C or lower

2.2 Reagents/supplies

- 2.2.1 *C. botulinum* Type C toxin, IRP 583
- 2.2.2 Pipettes, 1-mL, 2-mL, 5-mL, 10-mL, and 25-mL
- 2.2.3 Syringes, needle-locking, 3-mL, and 5-mL
- 2.2.4 Needles, 23-gauge x 1-inch
- 2.2.5 Serum vials, 50-mL, with caps and seals
- 2.2.6 M/15 phosphate buffered physiological saline with 0.2% gelatin

2.3 Test animals

- 2.3.1 Mink, young adult, from the same source and similar age. Vaccinate a group of 5 mink for each product tested. Three additional mink are held as unvaccinated controls.
- 2.3.2 Mice, 16-20 g, fifty mice are required to determine the LD₅₀ of the challenge inoculum. All mice must be from the same source colony.

UNCONTROLLED COPY

Supplemental Assay Method for Potency Testing of *Clostridium botulinum* Type C Bacterin-Toxoids

3. Preparation for the Test

3.1 Personnel qualifications/training

Technical personnel need to have a working knowledge of the use of general laboratory chemicals, equipment, and glassware; and have specific training and experience in the safe handling of clostridial toxins. Personnel need specific training in the care and handling of laboratory animals.

3.2 Selection and handling of test animals

3.2.1 Healthy, young adult mink are caged individually for the potency test.

3.2.2 Mink of either sex or any color may be selected.

3.2.3 Susceptible mink are required. They must not have been vaccinated with any *C. botulinum* Type C-containing product.

3.3 Preparation of supplies

3.3.1 Sterilize all glassware before use.

3.3.2 Use only sterile supplies (pipettes, syringes, needles, etc.).

3.3.3 Operate all equipment according to the manufacturers' instructions.

3.4 Preparation of reagents

3.4.1 M/15 phosphate buffered saline with 0.2% gelatin

| | |
|---------------------------------|-----------|
| NaCl | 8.5 g |
| Gelatin | 2.0 g |
| M/15 phosphate buffer at pH 7.4 | 1000.0 mL |

Stock solution of M/15 KH₂PO₄ (Solution A):

Dissolve 9.072 grams of KH₂PO₄ in one liter of distilled water.

Stock solution of M/15 Na₂HPO₄ (Solution B):

Dissolve 9.465 grams of Na₂HPO₄ in one liter of distilled water.

One hundred ninety-two (192) milliliters of Solution A are mixed with 808 mL of Solution B to provide 1000 mL of M/15 phosphate buffer at pH 7.4.

UNCONTROLLED COPY

Supplemental Assay Method for Potency Testing of *Clostridium botulinum* Type C Bacterin-Toxoids

Eight and one-half (8.5) grams of NaCl and 2 grams of gelatin are dissolved in the 1000 mL of M/15 phosphate buffer (pH 7.4). It is then autoclaved at $\geq 121^{\circ}\text{C}$ for 30 minutes. Store at 2° - 7°C for up to 6 months.

3.4.2 *C. botulinum* Type C toxin

C. botulinum Type C challenge toxin, IRP 583, is supplied by the Center for Veterinary Biologics (CVB). Each vial contains 1.3 mL. Vials may be combined to make sufficient quantity for each challenge. Store toxin vials at -70°C or colder.

4. Performance of the Test

4.1 Vaccination

Vaccinate the mink subcutaneously with the mink dose recommended on the label of the test serial final container using 5-mL sterile, disposable syringes with 23-gauge x 1-inch disposable needles.

4.2 Challenge of immunity

4.2.1 Challenge both vaccinated mink and unvaccinated control mink 21 to 28 days after the vaccination.

4.2.2 Dilute *C. botulinum* Type C toxin, IRP 583, 1:10 in M/15 phosphate buffered saline with 0.2% gelatin. Prepare the dilution by adding 1 mL IRP 583 to 9 mL diluent. IRP 583 diluted 1:10 contains $10^{4.0}$ mouse lethal dose fifty (LD_{50}) per 0.5 mL.

4.2.3 Inoculate the mink intraperitoneally with 0.5 mL of diluted IRP 583 using a 3-mL disposable syringe and 23-gauge x 1-inch needle.

4.2.4 Immediately following the mink inoculations, the 1:10 toxin dilution is further diluted to the final toxin dilutions of 1:100,000, 1:200,000, 1:400,000, and 1:800,000 in M/15 phosphate buffered physiological saline with 0.2% gelatin (see **Appendix**). Inoculate intraperitoneally each of 10 mice with 0.5 mL of each final toxin dilution. In most tests, at least 9 of 10 mice inoculated with the 1:100,000 dilution will die within 7 days and at least 9 of 10 mice inoculated with the 1:800,000 dilution will survive for 7 days.

Supplemental Assay Method for Potency Testing of *Clostridium botulinum* Type C Bacterin-Toxoids

5. Interpretation of the Test Results

Observe the mink for 7 days after challenge for signs of botulism and note deaths. For a valid test, the unvaccinated controls must die. If the test is valid and at least 80% of the vaccinated mink do not remain free of botulism, the serial is unsatisfactory.

Note: Moribund animals observed exhibiting clinical signs consistent with the expected disease pathogenesis that are unable to rise or move under their own power are humanely euthanized and considered as deaths as outlined in 9 CFR 117.4.

6. Report of Test Results

Report results of the test(s) as described by standard operating procedures.

7. References

Title 9, *Code of Federal Regulations*, part 113.110, U.S. Government Printing Office, Washington, DC.

8. Summary of Revisions

Version .03

- Section Leader and Director information have been updated.
- Toxin lot number and corresponding use dilutions have been updated throughout the document.

Version .02

This document was revised to clarify the practices currently in use at the Center for Veterinary Biologics and to provide additional detail. While no significant changes were made that impact the outcome of the test, the following changes were made to the document:

- The document number has been changed from BBSAM0213 to SAM 213.
- The Standard Requirement has been changed from 9 CFR 113.95 to 9 CFR 113.110.
- The Contact information has been updated.

UNCONTROLLED COPY


Supplemental Assay Method for Potency Testing of *Clostridium botulinum* Type C Bacterin-Toxoids

- **Section 4.3** information has been moved to **Section 5**.
- Toxin lot numbers and corresponding use dilutions have been updated throughout the document.

UNCONTROLLED COPY

Supplemental Assay Method for Potency Testing of *Clostridium botulinum* Type C Bacterin-Toxoids

Appendix: Mink Challenge and LD₅₀ Dilution Preparation

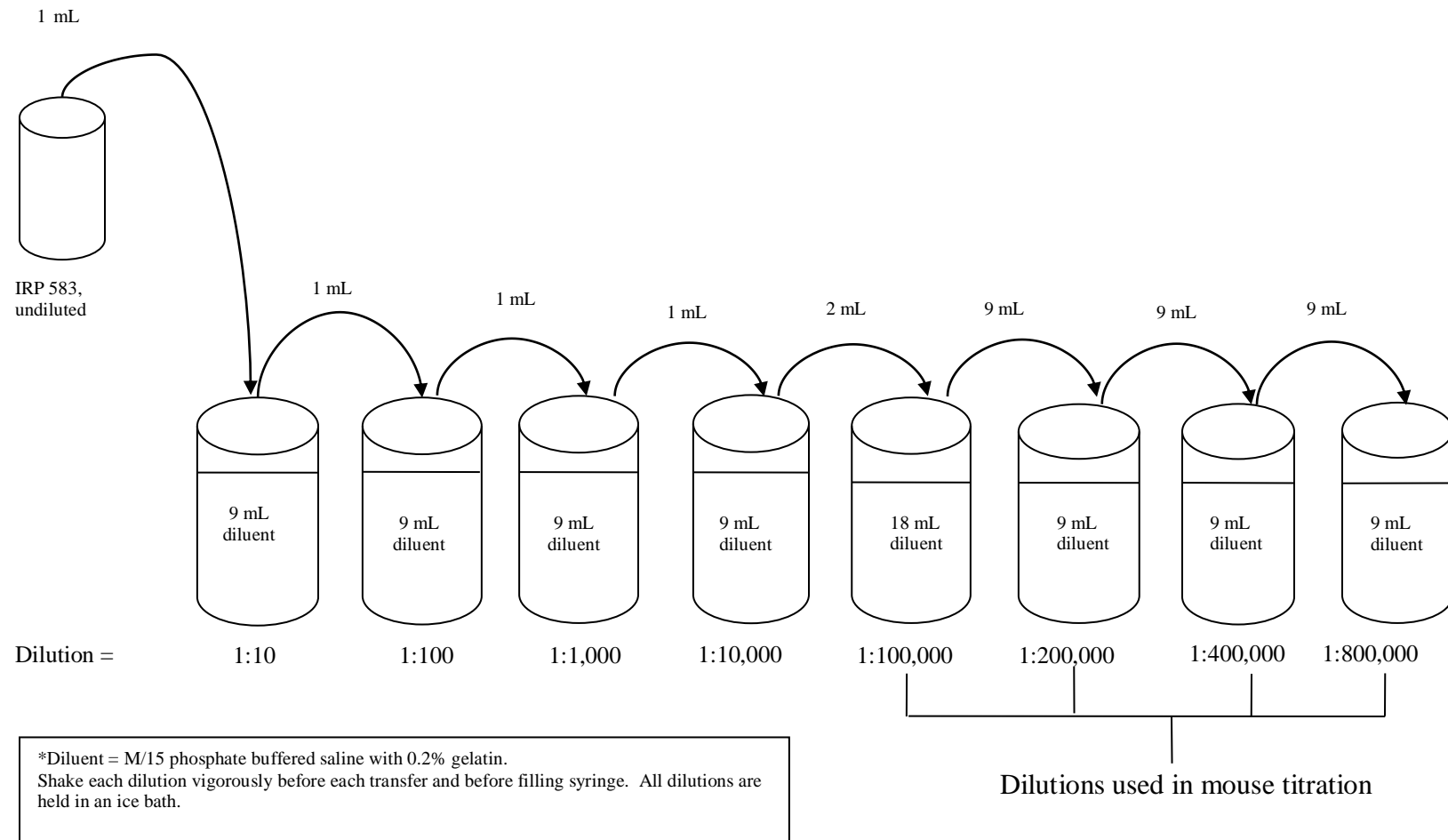

UNCONTROLLED COPY